



Charging KorePulse™

1. Connect the charging connector of the power adapter into the charging port of KorePulse™.
2. Plug in the power adapter to the desired outlet.
3. The LED indicators will blink red while KorePulse™ is charging.



* If KorePulse™ has been charged for three hours but not all the three battery LED indicators are on, unplug the power adapter to avoid danger and contact the customer service team.

* Do not use KorePulse™ while it is being charged.

4. When the three battery LED indicators have stopped blinking and shine continuously, KorePulse™ is fully charged.

* KorePulse™ typically takes 1.5 hours to full charge (when fully depleted).

5. When KorePulse™ has been fully charged, unplug the power adapter and remove the charging connector from KorePulse™.

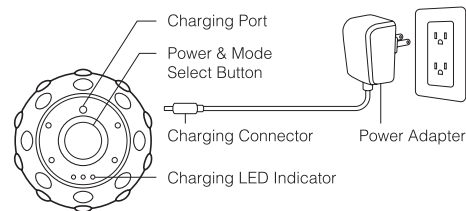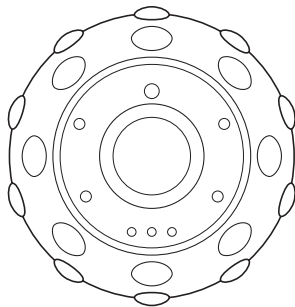



KorePulse™ has three different vibration speeds: Low, Medium, and High. It also features “Recovery Massage Mode,” which creates soothing pulses of vibration. Press the button to cycle between the different modes.

How to Use

- To turn on KorePulse™, press the button once. The first LED indicator light will turn on and the device will begin vibrating on Low.
- To increase the speed, press the button a second time. The second LED indicator light will turn on and the device will begin vibrating on Medium.
- To activate Recovery Massage Mode, press the button a third time. The third LED indicator light will turn on and the device will begin vibrating at variable levels.
- To increase the speed to the maximum level, press the button a fourth time. The fourth LED indicator light will turn on and the device will begin vibrating on High.

Turning Off KorePulse™

- Press the button a fifth time to turn off KorePulse™. The device will stop vibrating and all the LED indicator lights will turn off.
- To turn off KorePulse™ without cycling through all four vibration modes, press and hold the button for two seconds.

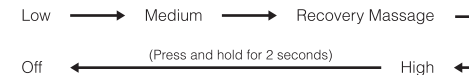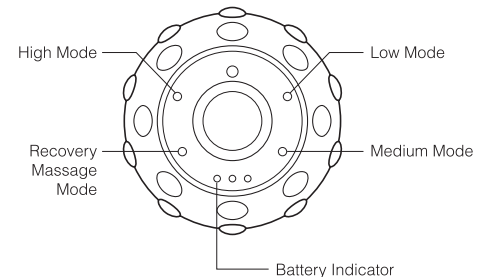

It's not recommended to use the device for prolonged periods of time. For safety, KorePulse™ will automatically shut off after 10 minutes.

Maintenance Instructions

- Before cleaning KorePulse™, make sure it is turned off and unplugged.
- Wipe KorePulse™ with a clean, soft cotton cloth.
- To avoid damage to KorePulse™, wipe the button area with a dry cloth.
- To avoid electric shock, keep KorePulse™ and its components away from liquids.

Proper Storage

- Keep KorePulse™ in a clean, dry place. We recommend putting it in its original packaging.
- Avoid exposure to direct sunlight or high temperatures.

Common Problems and Solutions:

Problem	Possible Cause and Solution
KorePulse™ does not vibrate when turned on.	Possible Cause The battery has been depleted.
	Solution Plug in KorePulse™ and charge it fully.
KorePulse™ vibration is too weak.	Possible Cause Low battery.
	Solution Please charge KorePulse™.
KorePulse™ stops vibrating just after a few minutes (even when fully-charged).	Possible Cause The surroundings are either too hot or too cold.
	Solution Only use KorePulse™ when the surrounding temperature is between 10°C and 40°C.
KorePulse™ suddenly stops working.	Possible Cause KorePulse™ is designed to automatically stop after ten minutes.
	Solution To continue using, restart the device.
KorePulse™ makes noise while working.	This is normal for the device.

Specifications:

Input	AC 100-240V, 50/60Hz, 0.5A MAX; DC 8.4V, 1000mAh
Battery spec.	DC 7.4V, 1200mAh
Charging time	About 1.5hrs (full charge)
Battery discharge time	50 mins
Working temperature	10°C-40°C
Maximum weight load	80kg
Product size	Ø95 x L90mm
Product weight	0.32kg
Material	Silicone and ABS
Product accessories	Power adapter, manual
Manufacturing origin	China

* The actual amount of time per battery charge depends on the number of uses, as well as the use environment.

- ☞ The device complies with the European Directives 2014/30/EU.
- ♻ At the end of its service life, this product may not be disposed of in the normal household waste, but rather must be disposed of at a collection point for recycling electrical and electronic devices. The materials are recyclable according to their labelling. You make an important contribution to protecting our environment by reusing, recycling, or utilizing old devices in other ways. Please ask the municipal administration where the appropriate disposal facility is located.
- Subject to change.

Precautions for use:

1. The device should only be used in a home or private environment.
2. For adult use only. Children should not play with the device.
3. This device may be used by persons with physical, sensory, or mental disabilities provided they are supervised.
4. Not suitable for pregnant women.
5. Do not let the charging port come in contact with moisture or submerge the device in water.
6. Do not modify or disassemble the device.
7. Do not use when the device is wet or damaged.
8. Consult a doctor before using this product if you:
 - Have a pacemaker;
 - Have any health problems;
 - Or are unsure whether this device is suitable for you.